EN

PL

BEAUTIFLY

B-Body

Blush

series

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


MEDICAL SILICONE
 PHOTON LIGHTS
 REDUCES CELLULITE
 LIGHT AND HANDY
 FIRMS THE SKIN
 IMPROVES THE ABSORPTION
 OF SLIMMING COSMETICS

The massage cushion is made of a skin-friendly, soft and safe medical silicone. A massage that stimulates circulation, exports toxins and other harmful substances of metabolism to remove unwanted cellulite and make the skin firmer. The technique of regulated vibrations in combination with photon light tightens the skin and allows it to maintain a healthy color. The product is ideal for use with a slimming cream or oils.

Technical data B—Body

Operating voltage:	3.7V
Product dimension:	110*110x78.8 (L*W*H)
Product weight:	~198g
Material:	ABS, medical silicone
Charging:	USB

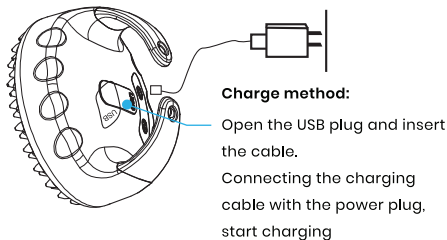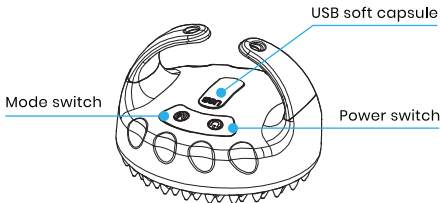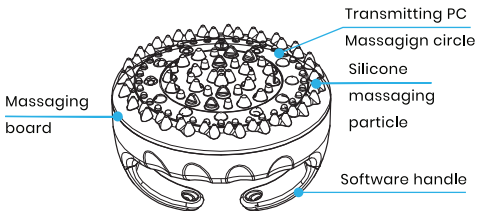Warning

- No use during charging.
- Do not disassemble or refit the product
- After usage, put in on the high place and keep it away from children.

Features

1. Firming massage.
2. Skin stimulation with special silicone insets.
3. Care and enhanced slimming effect thanks to photon lights.

Application method

Waist: Massaging upwards

Belly: Massagign in circularity

Hip: Massagign pulling

Thing: Massaging upwards

Arm: Massaging upward

Contraindications

1. Not suitable for pregnant or inflamed skin.
2. Keep the device away from children.
3. Do not use silicone or granular emulsion together with the product, so as not to damage the massage surface.
4. In case of any bad skin symptoms, consult your doctor.

How the B-Body device works

- 1 The device has 3 program
 - the first for red light,
 - the second for green light,
 - the third for blue light.
2. The device cannot be turned on while charging. When the power light flashes, it is charging. A steady power indicator light indicates charging is complete.

3. The device will turn off automatically after 15 minutes in operating mode. After being turned back on, it can be used continuously. Continuous blinking power indicator indicates low power requires charging.
4. The device has a memory function. After switching on, the device remains in the last used mode.
5. Switching on by long pressing the switch for 2 seconds. During operation in the M mode, 3 modes can be switched cyclically

Properties of photon lights

Red light reaches fibroblasts, which are responsible for the synthesis of collagen, elastin and hyaluronic acid.

It strengthens the activity of cells, speeds up metabolism, stimulates the skin to produce collagen proteins and delays the dying off of existing collagen. It also has a detoxifying and draining effect, increases the flow of the lymphatic system.

Blue light inhibits the development of inflammation.

It accelerates the healing of wounds, contusions and bruises and is used against actinic keratosis.

Green light acts at the level of the epidermis. Recommended for the reduction of discoloration and sun spots. Penetrates the basic layers of the epidermis, stimulates nutrients, helps regenerate mature skin.

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The crossed-out wheeled bin symbol indicates that in the European Union after the end of the product's useful life it must be disposed of in a separate, dedicated collection point. Do not dispose of these products together with unsorted general household waste

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English

Hereby, Beauty Factor Sp. z o. o. declares that the device equipment Beautifly B-Body is in compliance with Directive 2014/30/UE.

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TERMS OF WARRANTY

1. The product warranty is valid for a period of 24 months.
2. The warranty covers the period from the date of receipt of the device from Beauty Factor Sp. z o. o. and covers defects resulting from causes inherent in the received device.
3. The device will be repaired as soon as possible, not exceeding 21 working days from the date of acceptance of the device for repair by an Authorized Service Center.
4. The customer has the right to apply for replacement of equipment free from defects, if the Authorized Service finds that the removal of the defect is impossible or requires excessive costs.
5. If only part of the product is defective and can be detached from the part of the product working properly, the Customer's rights resulting from this warranty limit only to the defective part of the product.

In order to benefit from the warranty conditions, please send a complaint notification together with the warranty card to the e-mail address serwis@rgmedia.pl

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I represent that I have read and accepted the terms and conditions specified in the Warranty Certificate.