
GIO-450PINK

GIO-450BLUE

GIO-450YELLOW

PL

# innoGIO
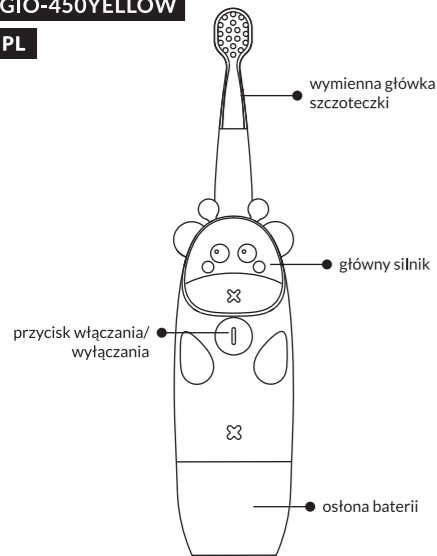
## GIO giraffe

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


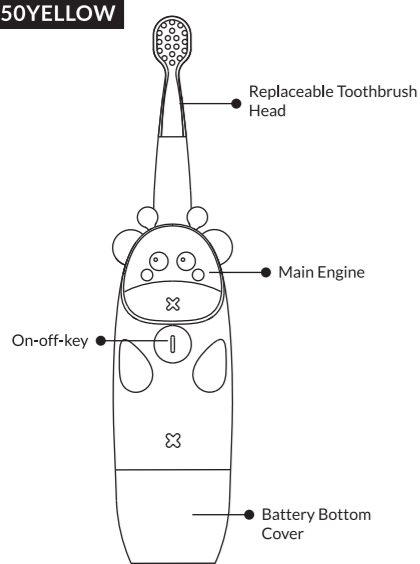
FIRST USE

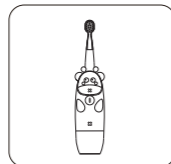
1. Take out the toothbrush.

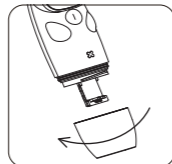2. Turn the bottom cover counterclockwise and open the bottom cover.

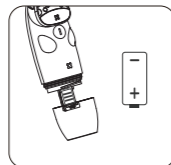
3. Take out the battery and the protective film, and replace the battery (Positive down).

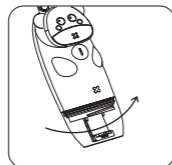
4. Put back the bottom cover and screw it clockwise.


USAGE

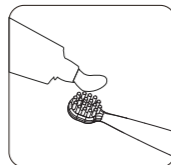
1. Squeeze toothpaste on the brush.

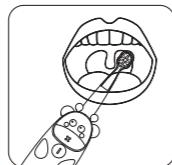
2. Put the toothbrush in your mouth.

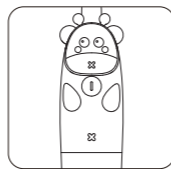3. To select an appropriate gear, gently press the switch. The timer counts down 2 minutes of brushing. There is also a reminder for changing the brushing area every 30 seconds. Press and hold for 3 sec. to turn it off.

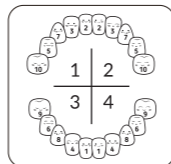4. Brush in line with dentists' recommendations: each tooth area shall be brushed for 30 sec. for cleaner teeth.

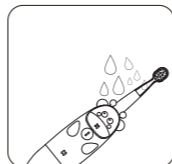5. Rinse your mouth with water and clean the toothbrush with water after brushing your teeth.

PRODUCT INFORMATION

Name: InnoGIO GIOgiraffe Sonic Toothbrush For Children  
Model: GIO-450PINK / GIO-450BLUE / GIO-450YELLOW  
Power supply: 1x AAA battery (no battery included)  
Product size: 160.5 x 35 x 35 mm  
Rated voltage: 1.5 V  
Rated power: 0.51 W  
Rated current: 0.34 A  
Waterproof grade: IPX7

FUNCTIONS

FOUR GEARS

- First Gear Seven-color Mode: turn the device on and it will display different colors: red, orange, yellow, green, blue, indigo and purple, the interval of color change is 0.5 sec. without vibration, the total brushing time is 1 minute. This gear is recommended for children who do not tolerate vibrations.
- Second Gear Gentle Brushing Mode: Second Gear Gentle brushing Mode with 12,000 pulsating strokes per minute: the light remains unchanged for 2 minutes. This gear is recommended for children aged 2-5 years.
- Third Gear Standard Brushing Mode: Third Gear Standard Brushing Mode with 16,000 pulsating strokes per minute: the light changes to the next color after the second gear, the total brushing time is 2 minutes. This gear is recommended for children aged 5-8 years.
- Fourth Gear Strong Brushing Mode: Fourth Gear Strong Brushing Mode with 20,000 pulsating strokes per minute: the light changes to the next color after the third gear, the total brushing time is 2 minutes. This gear is recommended for children aged 8-12 years.

2-MINUTE TIMING

For the second, third and fourth gear modes, after 30 sec., 60 sec. and 90 sec., the lights flash once, twice and three times respectively, and the engine stops for one time, two times and three times, respectively. The flashing and stopping time is 0.5 sec. to remind the user to change the brushing area 120 sec., the light flashes twice, and then it will turn off.

MEMORY FUNCTION

The device can enter the memory mode in every gear. Press and hold to turn off or shut down automatically. Next time you turn the device on, it will automatically enter the last used.

CHARACTERISTICS

**DuPont Brush**  
The brush is made of DuPont material, with a rounded end, which is strong and wear-resistant. It protects the baby's gums while cleaning teeth efficiently.

**Colorful Lights**  
The toothbrush has a cute animal design with a colorful glare - children will love it, while learning the good habit of teeth brushing.

**Silicone Sleeve Design**  
The unique and creative silicone design is three-dimensional, providing a more comfortable and smoother feeling. Kids will love the feeling of tooth brushing.

**IPX7 Waterproof**  
The waterproofing grade of the device is IPX7, which means the toothbrush can be washed entirely, and even used during bathing.

CONTENTS OF THE PACKAGE

- Toothbrush
- Replaceable brush head 2 pcs
- Instruction

ATTENTIONS

1. The normal working temperature is 0-50°C.
2. If the device is not used for a long time, please remove the battery, and store the toothbrush in the box, in a cool place.
3. Leave the toothbrush head to dry after brushing teeth.
4. When the product is in use, it reaches its maximum output power. The product mustn't work continuously for a long time. Generally, each use time should not exceed 10 minutes.
5. Never walk or run when using the toothbrush.
6. Make sure that children are always supervised by an adult while brushing their teeth.

innoGIO

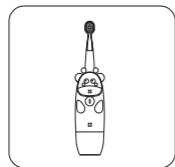
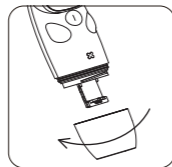
Produced in China for Innogio sp. z o.o.  
Ostródzka street 74H, 03-289 Warsaw, Poland  
www.innogio.com  
Date of issue: 05/2021

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