TAKE CARE!

babuŇo

<b>GB</b>	Electronic scale	2
<b>PL</b>	Waga elektroniczna	3
<b>BG</b>	Електронна везна	4
<b>CZ</b>	Digitální váha	5
<b>DE</b>	Elektronische Wiege	6
<b>ES</b>	Báscula electrónica	8
<b>EE</b>	Elektroniline imikukaal	9
<b>FR</b>	Balance électronique	10
<b>GR</b>	Ηλεκτρονική ζυγαριά	11
<b>HR</b>	Elektronska vaga	13
<b>HU</b>	Digitális babamérleg	14
<b>IT</b>	Bilancia elettronica	15
<b>LT</b>	Elektroninės svarstyklės	16
<b>LV</b>	Elektroniskie svāri	17
<b>NL</b>	Elektronische weegschaal	18
<b>PT</b>	Balança eletrónica	20
<b>RO</b>	Cântar electric	21
<b>RU</b>	Весы электронные	22
<b>SE</b>	Elektronisk babyvåg	24
<b>SK</b>	Digitálna váha	25
<b>UA</b>	Весы электронные	26

## Electronic scale

Read the user's manual and warnings before using the device for the first time. Pay special attention to safety tips to avoid accidents and/or damage to the product when using it. Keep the user's manual for reference purposes!

The modern electronic scales are a device for checking the baby's weight. The device is especially useful for weighing infants and premature babies, where the weight is an exceptionally important development factor.


With your needs in mind, we created a product that is easy to use, accurate and highly reliable. The scales have a taring function across the entire weighing range. The scales automatically lock the result for easier reading. We used our best efforts to make sure the product is of top quality, comfortable and safe during weighing. The weight limit (max. 30 kg) is designed for babies.

This electronic scale is a device for checking the baby's weight at home, it is especially useful for weighing infants and premature babies, where the weight is an exceptionally important development factor. We do our best efforts to ensure the product with top quality, comfortable and safe during weighing.

### TECHNICAL DATA

- Device dimensions: 700 x 319 x 42.5 mm
- The Capacity: 40 g ~ 30 KG
- Press on and auto off
- Division: 2G
- 4 high precision strain gauge sensors
- Tare function
- Low battery and over load indicate
- Auto lock/hold function (easier weighing of moving babies, the result lock on display automatically within 1-2 seconds)
- Large LCD: 81 x 16.9 mm with backlight
- Battery: 3 x 1.5 V AAA batteries

### DEVICE INSTRUCTION
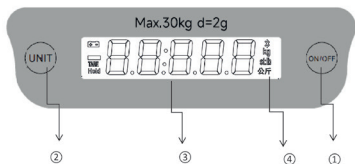
#### 1. ON/OFF button - tare button

- Press the button to switch on or off the device
- While use a towel for more comfortable weighing, first put the towel on the platform and then press the button to set the display to "0.00".

#### 2. Unit button

- Press the button to convert the units of measurement: gram (g), kilogram (kg), pound (lb) or stone (st:lb).

3. LCD screen - Display weight.

#### 4. G/KG/LB/LB:ST indication

-The unit digits comes on when gram (g), kilogram (kg), pound (lb) or stone (lb:st) is set as unit of measurement.

### BATTERY INSTALLATION

1. Before using the scale, check the battery compartment and remove any plastic insulation strips
2. Insert the batteries (3 x 1.5V AAA) into the battery compartment
3. If you are replacing the batteries, make sure the positive and negative contacts are properly aligned.

### DEVICE OPERATION

#### WEIGHT

1. Please use the scale on the hard & flat surface. Do not weight on carpet or any soft surface as it may affect accuracy.
2. Press "ON/OFF" button, the display shows "Hi";
3. After display "0.00 kg", then ready for weighing;
4. Put baby or other stuff on the scale platform, and LCD display will show the weight;
5. Press "ON/OFF" key to turn off the scale, or scale will be auto off after 30 seconds while in idle.

#### AUTO LOCK FUNCTION

When scale has the weight of the baby, do not move the baby off the platform, the result lock on display automatically within 1-2 seconds.

#### TARE FUNCTION

It could weight continuously:

- First, put the first baby on the scale platform to get weight, take away the baby, the baby's weight still on the LCD.
- Second, press "Tare" button, LCD will turn to zero, then put the next one on the scale for weighing. ( this function is repeatable)

#### UNIT CONVERSION

Press "unit" button to convert between "g"/"kg"/"lb" and "st:lb".

#### OVERLOAD INDICATION

When the weight is over the "max capacity" LCD display will show "0-Ld"

#### LOW POWER INDICATION

When LCD display shows "LO", Please replace the new battery with correct polarity.

#### SAFETY PRECAUTIONS

- Ensure to place the scale on flat, dry and hard surface when weighing.
- Use this device for its intended purpose at home.
- Do not load object that over the scale's max capacity: 20kgs.
- Avoid shaking or hitting the scale.
- Do not submerge it into water.
- When in storage, do not place anything on the scale.
- The optimal temperature for operating and storing the scale is from 5 to 30 degree celcius.

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